



The BPL Video Laryngoscope

A new age Video Laryngoscope for a
fast and effective intubation



Follow Us on



www.bplmedicaltechnologies.com



Happier Living Everyday

The BPL Video Laryngoscope

Compact, Lightweight & Portable

4" INCH LCD DISPLAY

High resolution 4' anti-Reflective LCD display, designed for maximum visibility

180° SCREEN ROTATION

180° Horizontal Screen Rotation along with Vertical movements for enhanced clinical assistance


CMOS CAMERA WITH LED

2MP CMOS camera, High Definition Anti Fog lens and integrated LED light source for high performance visualisation



QUICK DATA TRANSFER

Videos and images can be transferred via a micro USB port for easy storage and access

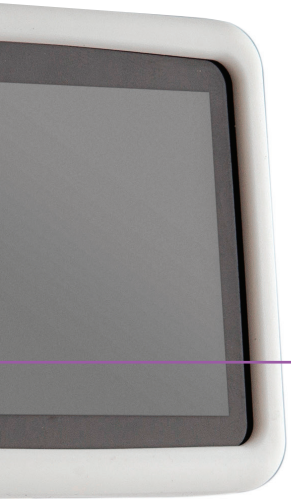

BATTERY INDICATOR

Rechargeable Lithium ion battery with a 2Hr capacity with a low battery indicator



BUILT IN MEMORY

4GB built in memory, for continuous video recording and image capturing with an use of a button.

ERGONOMIC HANDLE

Anti microbial, ergonomically designed, reusable handle to assist clinicians while use



DISPOSABLE BLADES

Affordable, sterile, vacuum packed disposable blades connects securely to the video laryngoscope handle with a lock. Available in 3 Sizes MAC2, MAC 3, MAC4



Product Specifications

Display	Size Monitor Resolution Ratio Aspect Ratio Video Refresh Rate	4" Full View LCD 800x480 (RGB) 4:3 30 FPS
Blades	Disposable Blades Available Sizes	Mac2, Mac3, Mac4
Laryngoscope Blade Camera	CMOS 2Mega Pixel High Resolution Resolution ratio Field Angle Illuminance	2.0M Pixels 66° ≥800 LUX
Handle	Anti-microbial, Comfortable ergonomic design	
Transport/ Storage Condition	Temperature Humidity Atmospheric Pressure	-10 to 45°C ≤93% 500hPa to 1060hPa
Working Condition	Temperature Humidity Atmospheric Pressure	-10 to 45°C 30% - 85% 700hPa to 1060hPa
Battery	Type Working Time Voltage Capacity Battery Life Cycle Charging time	Rechargeable Lithium ≥ 130 mins 3.7V 3200mAh ≥ 300 times 4-6 hrs
Power	Input Output	AC 100-250V, 50Hz DC 5V, 2000mA
Order information	BPL Video Laryngoscope kit Disposable Blade MAC2 Disposable Blade MAC3 Disposable Blade MAC4	91MED170 (includes 4 disposable blades) 91MED175 91MED176 91MED177

CERTIFIED ISO 13485 : 2016 COMPANY

BPL Medical Technologies Private Limited

Regd. Office: 11th KM, Bannerghatta Road,
Arakere, Bangalore - 560076, India.

Toll Free: 1800-4252355

Website: www.bplmedicaltechnologies.com

For Enquiries: sales.medical@bpl.in

CIN: U33110KA2012PTC067282

