

Respirare Technologies

Respirare VT 200 BiPAP device gives you the flexibility to adjust the settings and accommodate various patients conditions.



Scan Report

Simplest way to get one day compliance report by scanning the QR code.



In.trigger

In.trigger lets you set min and max limits on either side of the patients ideal inspiratory time.

Technical Specification

Modes	CPAP , Auto CPAP, T, S, ST, PC, VAPS mode with T, ST, PC
Operating Range	IPAP: 6-30 cmH2O EPAP: 4-20 cmH2O in T, S, ST, PC mode 4-20 CmH ² O in CPAP mode VT:120-1200 ml in VAPS mode
Sound Pressure Level	<28 dBA with uncertainty of 2 dBA
Leak Compensation	10-40 LPM
Peak Flow	200 LPM
Respiratory Rate	4-60 BPM
Ramp/Rise Time	Auto, 0-60 mins/ 1-5 Level
In.trigger	1-10 level
Humidification Level	Off, 1-5 level
Alarms	High Pressure, Low MV, Low TV, Mask & Leak Alert
Data Management	Yes (1 day, 7 day, 31 days) of detailed data stored in the device Scan to see one day compliance report
Dimensions (LxWxH)	290x160x110 mm (Device with Humidifier)
Weight	1800 g (Device with Humidifier)
Power Adapter	Input Range : 150-240V, 50/60 H ² , 5A 60 Watt (typical Power Consumption)
Operating Temperature/ Humidity/ Altitude	+5 c to +35 c/ 10-95% relative humidity, non-condensing/sea level to 2,591 m, air pressure range 1013 hpa to 735 hpa
Storage and Transport Temperature/Humidity	-20 c to +60 c / 5-95% relative humidity, non-condensing
Supplement Oxygen	Maximum Flow: 15L/min (S, ST, T, PC) 4L/min (CPAP)
Standard Air Filter	Pollen Filter (Washable)
Tubing	Standard (19 mm)
Air Outlet	22 mm Air Outlet (STD)
Warranty	2 Years (Replacement)*
IEC	IEC 60601-1-1 Class 2 (double insulated, Type BF, Ingress Protection IP22) IEC 60601-1-2 IEC 80601-2-12

*Term & Conditions Apply

Key Features

VAPS (Volume-Assured Pressure Support)

VAPS mode mitigates hypoventilation by providing a constant tidal volume, the volume of air inspired or expired in a single breath during regular breathing. It automatically adjusts pressure support to maintain an optimal consistent tidal volume in which IPAP automatically increases or decreases to maintain an idealized target tidal volume.

i-Mode

This feature will turn the device **On** for you have the mask on without the need to press any buttons and will also turn the device **Off** once you take your mask Off.

A-Flex

The A-Flex is a comfort feature that will lower the pressure when you exhale and slowly increase the pressure when you inhale and will be available in Auto CPAP mode only.

Product Codes

Respirare BiPAP with Humidifier

VT 200

Note: Technical specification may change without notice. Always refer to your device user manual.

Manufactured By:
Deck Mount Electronics Pvt. Ltd 260, Phase IV,
Udyog Vihar, Sector 18, Gurugram, Haryana 122015

Saving millions of lives


Deck Mount BiPAP offers the best solution for patients who need extremely reliable, non-invasive type of ventilation either at home or in hospital care, and is equipped with all the latest innovative technology features. Deck Mount BPAPs are an alternative when a CPAP machine isn't working or well tolerated. They deliver a higher air pressure on inspiration, and drop the pressure on exhalation. This makes them more effective or more comfortable for some people with obstructive sleep apnea.



Hospital



Home

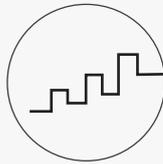

Ramp



In.trigger



Humidity



2 Years Warranty



wecare@deckmount.in



www.deckmount.in



Toll Free No.
1800 30 92 499



Deck Mount Electronics
Pvt. Ltd, 260, Udyog Vihar,
Phase 4, Gurugram 122015,
Haryana, India



Make In India

